


# **MASSAGE CHAIR**

## USER MANUAL

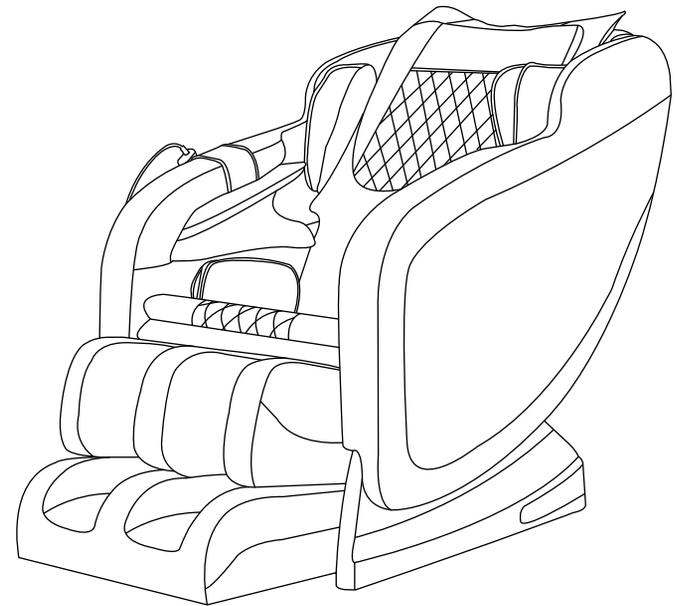

**Favor-MM550**

PLEASE READ AND SAVE THESE INSTRUCTIONS

Thank you very much for purchasing our massage chair. This operating manual contains important information about how to use the chair safely and correctly. Be sure that you have read all safety precautions thoroughly before using your chair. Keep this manual nearby so that you can refer to it anytime.

We reserves the right to change the product design at any time without notice. Pictures are for reference only and may vary from actual product.

## CONTENTS

Safety Precautions and Maintenance

Components and Description

Main Features

Operations

Specifications

## NOTICE

Our massage chair is 100% tested before leaving the factory.


## IMPORTANT

1. Make sure the chair is connected to an outlet having the correct voltage.  
Make sure that the power cord is seated fully into the power outlet  
Otherwise it may result in fire or electric shock.
2. Do not modify the plug provided with the product if it will not fit the outlet.
3. Do not use the chair outdoors.
4. Do not use the chair in any way not prescribed by the Operating Manual.
5. Do not install the chair in an area of high humidity such as a sauna.
6. Persons with the following conditions should seek medical advice before using the chair:
  - Those in the early stages of pregnancy or immediately after giving birth.
  - People with weak bones, osteoporosis, or spinal fractures.
  - People with sensory impairment caused by circulatory or other serious ailments.


- People suffering from heart disease.
  - People have a pacemaker or other medical devices prone to electronic interference.
  - People who are heat sensitive.
  - People under a doctor's care, or those experiencing unusual physical discomfort
  - Persons who are suffering from any other serious impairment or illness.
7. Keep massage sessions to no more than 30 minutes and rest after every massage.
  8. If you begin to experience pain or other physical discomfort while using the chair, immediately discontinue use .
  9. Do not fall asleep in the chair during a massage.
  10. Do not use the chair after consuming alcohol.
  11. If the chair does not start, or if there is any type of malfunction, please stop the chair immediately and turn the main POWER switch to the OFF position.
  12. Contact us for assistance regarding professional repairs or maintenance.
  13. Make sure to turn the main POWER switch to OFF after every use and to unplug the chair from the power outlet.
  14. Make sure to turn the main POWER switch to OFF and unplug the chair from the power outlet during cleaning.
  15. Never plug or unplug the chair from a power source (surge protector, wall outlet, ect.) with wet hands.
  16. Never use a wet cloth to clean the power switch, power cord, or any other electronic component.

## COMPONENTS AND DESCRIPTION

Chair:

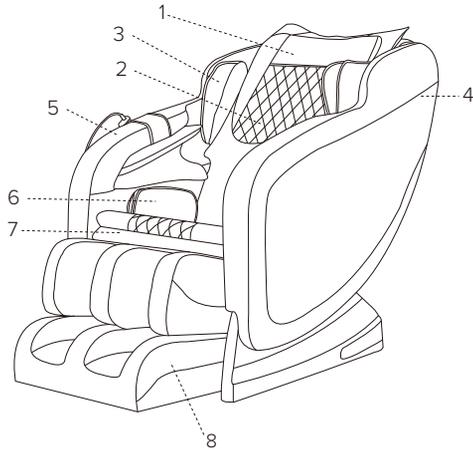


- 1.Pillow
- 2.Backpad
- 3.Shoulder Units
- 4.Backrest
- 5.Armrest
- 6.Hip Airbags
- 7.Seat
- 8.Footrest

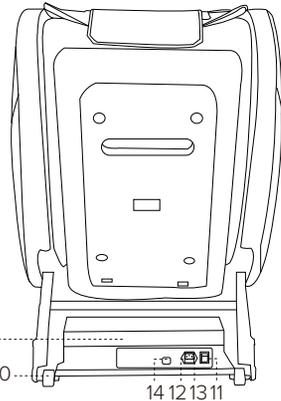

- 9.Power Supply
- 10.Wheels
- 11.Power Switch
- 12.Electrical Fuse
- 13.Cord Connector
- 14.Remote Control

## REMOTE CONTROL

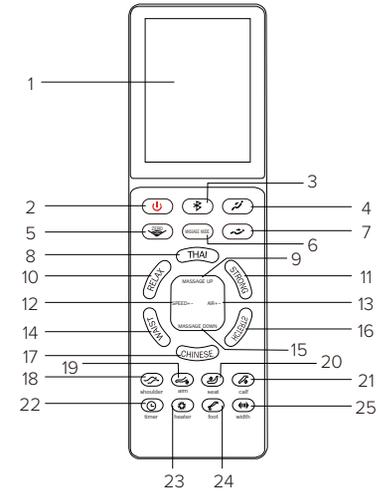

1. TFT Display
2. Power Button
3. Bluetooth On/Off
4. Backrest Rise
5. Zero Gravity
6. Massage Mode
7. Backrest Down
8. Thai Massage Mode
9. Massage Roller Up
10. Relax Massage Mode
11. Strong Massage Mode
12. Massage Roller Speed +-
13. Airbag Pressure +-
14. Waist Air Pressure On/Off
15. Massage Down
16. Stretch Massage Mode
17. Chinese Massage Mode
18. Shoulder Air Pressure On/Off
19. Arm Air Pressure On/Off
20. Seat Air Pressure On/Off
21. Calf Air Pressure On/Off
22. Timer On/Off
23. Heater On/Off
24. Foot Roller On/Off
25. Massage Roller Width +-

## ACCESSORIES FROM BOX



## INSTALLATION INSTRUCTIONS



First, lift the backrest.



When you hear a click, the backrest is installed.



Install the backpad



Connect backpad zipper



## INSTALLATION INSTRUCTIONS



Fixed pillow



Connect the heating power cord



Connecting the trachea


## INSTALLATION INSTRUCTIONS




Connect the hook to the host



Install arm screw

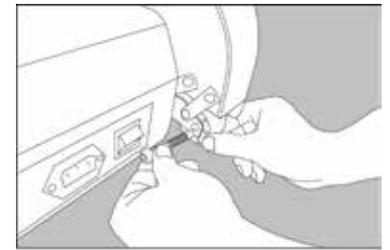

Fix the wheel with a screw.

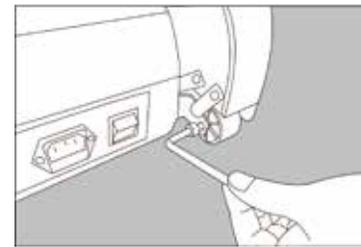


Fix the wheel with a screw.



Plug in the remote control wire.

## AUTO MODE

### Select Auto massage mode

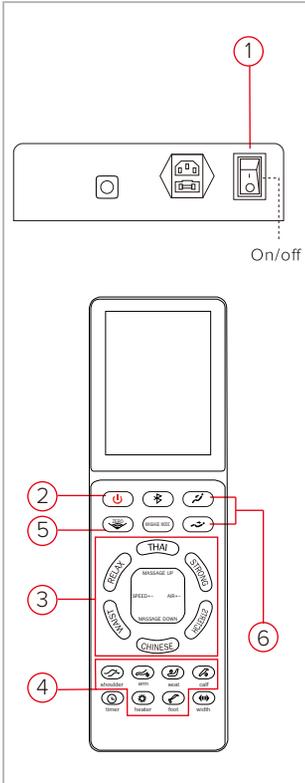

1. turn on the power:when the power is turned on, press UP and DOWN to lower or raise backrest.
2. press the power button:when press power button, it will automatically start Auto program one. turn off :press and hold POWER button for 2seconds,the unit will turn off.
3. you can switch auto program by press auto program buttons,we have 6 automatic modes:

- relax
- thai
- strong
- stretch
- chinese
- sleep



- you can switch to next auto program anytime, don't have to wait until current program finish.
4. in any auto program,you still can turn on/off: Shoulder airbags,arm airbags,seat airbags, calf airbags,heater and foot roller.
  5. Zero:In any auto programs,press ZERO into zero gravity position.
  6. UP and DOWN:In any auto programs,press UP and DOWN to lower or raise backrest.
  7. timer:press POWER button,you can switch time from 5mins-30mins.

## CUSTOM MESSAGE MODE

- By press "MESSAGE MODE", auto massage will stop  
You can turn on a single massage area, or multiple areas.
- Switch the massage mode by press "MESSAGE MODE"
- Adjust the width of back rollers by press "Width"
- Adjust the position of back rollers by press "MESSAGE UP and MESSAGE DOWN"

### About the remote control panel

- BLUETEETH:Press to turn on or off Bluteeth
- UP and DOWN:Press to lower or raise backrest and feet
- ZERO:Press to enter zero gravity mode
- MESSAGE MODE:Press to enter custom massage mode
- AUTO1-6:press to enter differents auto programs
- MESSAGE UP and MESSAGE DOWN:press to adjust the back roller position(In custom mode)
- SPEED+/-:Press to adjust the speed of back rollers
- AIR+/-:Press to adjust the air pressure of the airbags
- Shoulder:Press to turn on/off Shoulder airbags
- Arm:Press to turn on/off Arm airbags
- Seat:Press to turn on/off Seat airbags
- Calf:Press to turn on/off Calf airbags
- Timer:Press to turn on/off timer
- Heater:Press to turn on/off the heater
- Foot:Press to turn on/off foot rollers
- Width:Press to adjust the width of back rollers(In custom mode)



### Leg extension settings


- Footrest unit will automatically extend to suit taller individuals.

## AFTER USE

### RESET THE CHAIR



Make sure children, pets, and foreign objects are away from the chair.

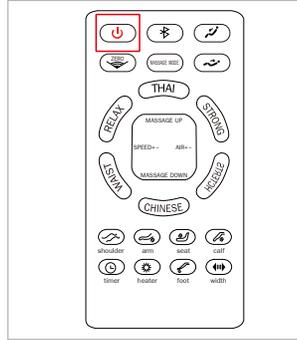

Press and hold power button for 2seconds, the reclining seat will return to its storage position and power off.


## MAINTENANCE

1. Anyone other than an authorized service person must not repair or disassemble the chair.
2. Stop the main power after massage is finished.
3. If you are not going to use the product for some time, store the power cord in a dry and clean place.
4. Do not leave the chair outdoors, in direct sunlight, in a high temperature environment, or near fire.
5. Use a piece of soft, dry cloth to clean the chair. Do not use alcohol or household cleaners on the chair.
6. Do not operate the chair continuously for long periods of time.

## TROUBLE SHOOTING

1. Check the power connector is fully plugged into the wall outlet if the chair doesn't work after the main power is connected.
2. The sound of the motor during operation is a normal part of the product's operation.
3. The chair will automatically power off after a massage. Press the power button to return the chair to an ON state.

### THE MAIN POWER

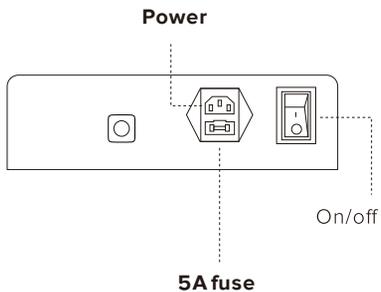

Shut down the main power

## Favor-MM550

### FUNCTIONS

---



Auto



Arm



Foot



Air Squeeze



Zero-g



Heating

---

Measurement: 46.5" \* 25.6" \* 29.5"    52.8" \* 25.2" \* 9.8"

G.W/N.W:149.8LB/126.7LB    42.1LB/33.7LB

Rated Power:100W

Rated Voltage:110V-60Hz

